



CONFIANZA[®] PRO 8-20

PCI GUIDE WIRE



YOUR DREAMS. WOVEN TOGETHER.

CONFIANZA PRO 8-20

— High penetrability and straightness —




0.008" Tapered Tip

Improved penetration force



20g Tip Load

Helping to cross the fibrous cap and the calcified lesion

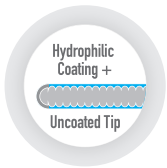


Uncoated Tip & Hydrophilic Coating

Catching the entry point and enhancing the penetration efficacy




1mm 45° Pre-shape

Less likely to lose the shape during the procedure

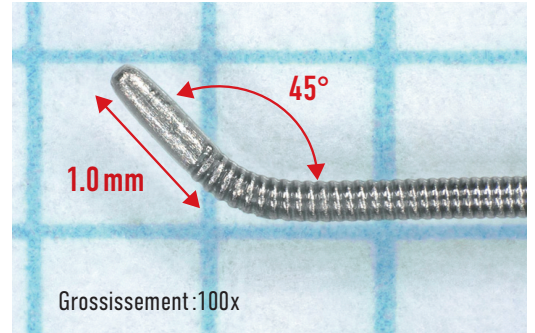


Tip load (gf)

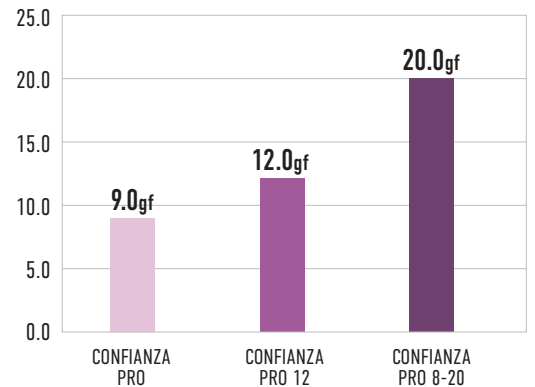

Structure & Ordering Information

Product	Catalog No.	Diameter	Usable Length	Spring Coil Length	Radiopaque Length	Tip Shape
	AGH143192	0.36 / 0.20mm (0.014 / 0.008inch)	190cm	17cm	17cm	Straight
	AGH143192P					Pre-shape
	AGH143392		300cm			Straight
	AGH143392P					Pre-shape

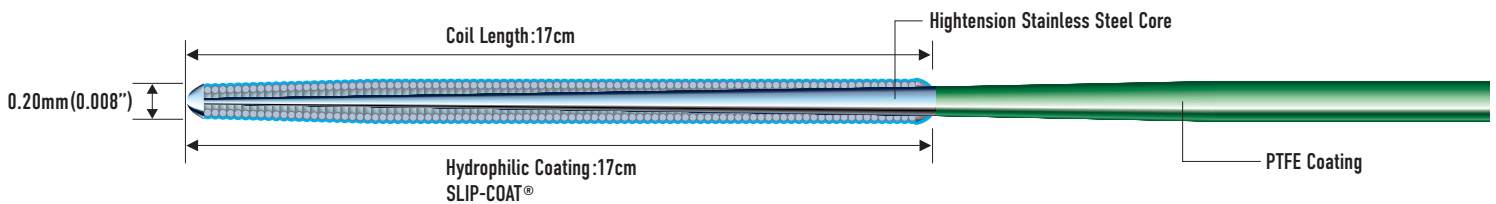

Toll Free: 855-286-9473 Fax: 949-756-8146
WWW.ASAHI-INTECC-US.COM



for Android



for iOS

ASAHI product information

© 2022 ASAHI INTECC CO., LTD.
"CONFIANZA" is trademarks or registered trademarks of
ASAHI INTECC CO., LTD. in Japan and other countries.
*Coated with SLIP-COAT® coating.
K22242_B_LC_1.0
MKT0827 Rev.A